

Mobile apps to record Sleep Data

Mobile Apps

Search for **SleepImage** and add






Note: -

- SleepImage mobile App is downloadable from either the Google Play Store for android phone (To enable location) or the Apple App Store for IOS phone.
- Recorder is paired with the App using Bluetooth

Pairing the Recorder with the App

Need to pair the recorder (Checkme 02) with mobile phone for data to be recorded:


- 1 Turn on mobile device's Bluetooth & location
- 2 Open the Sleep Image App
- 3 Tap on Setting 
- 4 Connect New Recorder
- 5 Press search
 - o Mobile phone will search for the device
- 6 Identify device serial number (last 5 digit)
- 7 Tap Setting 
 - o Sign in
- 8 To register patient to the mobile phone recorder for recording sleep data

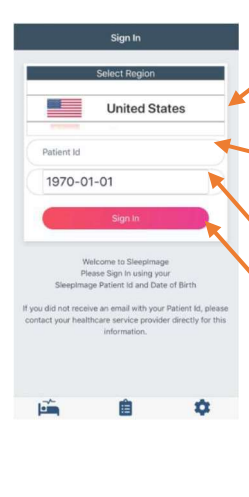

The screenshot shows the 'Sign In' screen of the Sleep Image app. It features a 'Select Region' dropdown menu with 'United States' selected. Below this is a 'Patient Id' input field containing '1970-01-01'. A pink 'Sign In' button is positioned below the input fields. At the bottom, there is a navigation bar with icons for home, list, and settings. To the right of the screenshot, four orange arrows point to the 'Select Region' dropdown, the 'Patient Id' field, the 'Sign In' button, and the 'Sign In' text in the instructions below.


- Choose Geographical Region Other
- Enter Patient ID from SleepImage Website generated
- Enter Date of Birth
- Sign In
If patient ID & DOB does not meet, it won't register

- 9 Successful registered

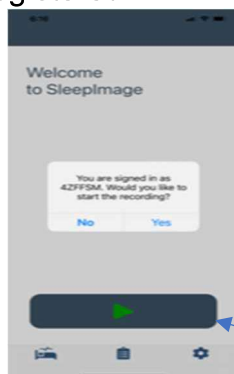

Press Start

10 Mobile phone will show the reading same as recorder.

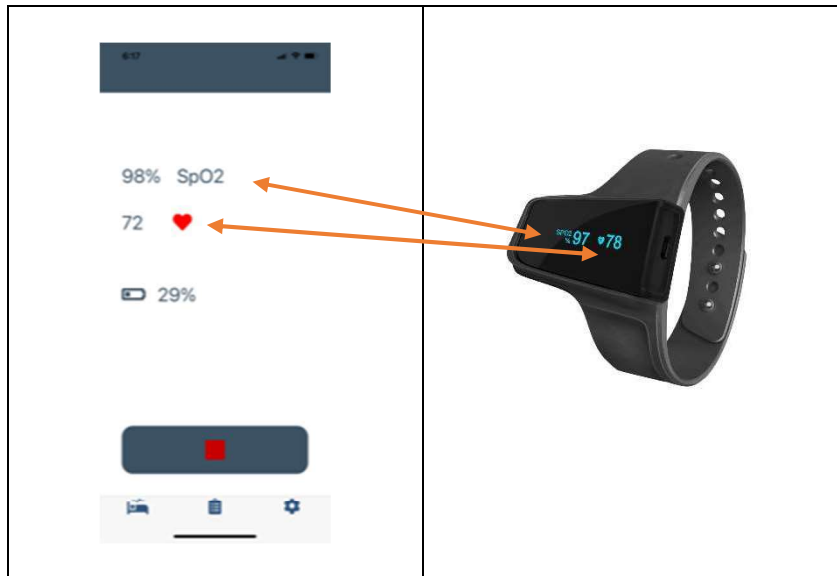

Note:

- Both reading for SpO2 and heart rate must be the same.

11 Press stop next morning when wake up